




HAELO

Quick Start Guide

HAELO
symphony one™

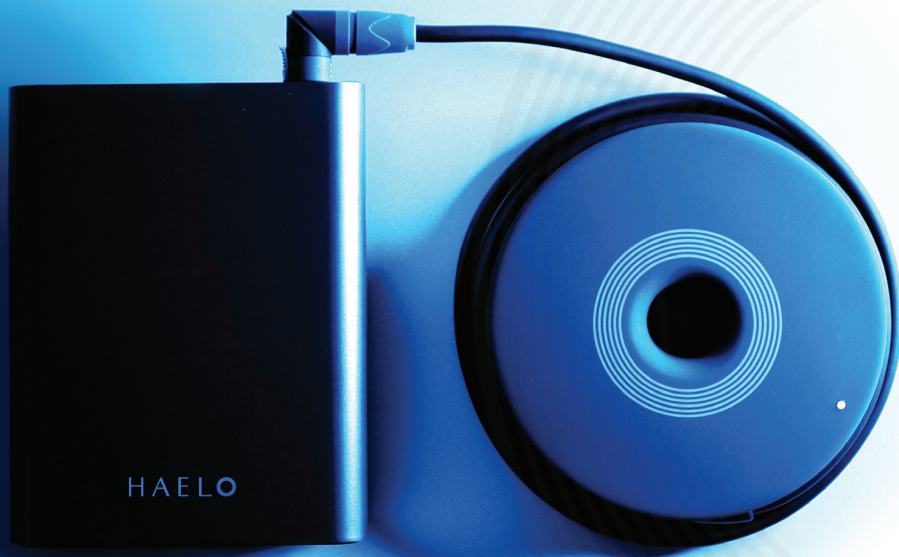



The Future of Symphonic Frequency has Landed!

Finally, a sound frequency therapy device that truly represents next-gen PEMF technology delivered through a patent-pending system that exudes simplicity, supports peak performance and an optimized life!

HAELO was born out of the desire to meet an active market that yearns for easy to use, plug and play technology that also delivers on wellness goals beyond expectation.


World class athletes are benefiting from HAELO's proprietary technology that delivers symphonic sine waves through a magnetic coil creating an entirely new realm of possibility for peak performance and faster recovery.

Designed with the active person in mind, HAELO offers thoughtfully arranged frequency sets that deliver what you need, when you need it—whether at home or on-the-go.

No matter what your goals are, HAELO will meet you with its portable travel-friendly design and smart app giving you the tools you need to go the extra mile to excel.

Energize your body, sharpen your mind, fuel your life! Tell us, how do you HAELO?


Getting Started

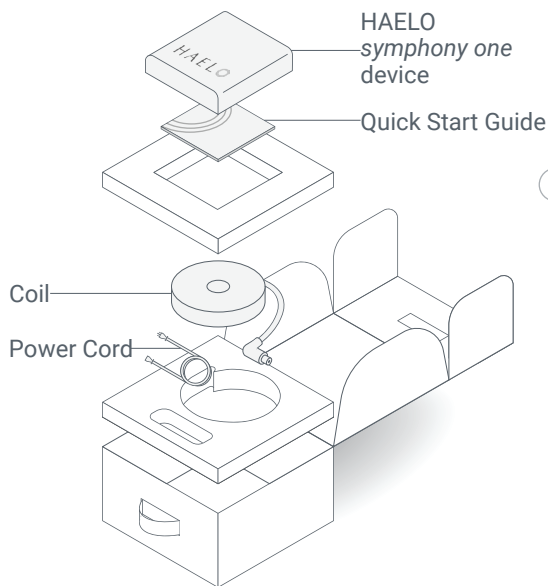


- 1 Remove HAELO from the protective foam.

Lift the protective foam layer to remove the coil and power cord.

Keep your HAELO shipping box as it's designed to travel and serves as a protective case.

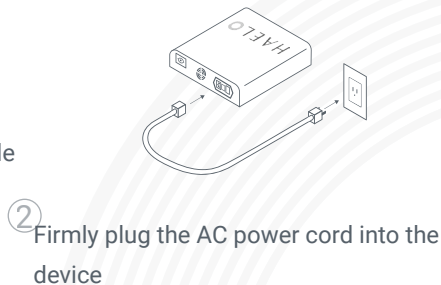


- 2 Firmly plug the AC power cord into the device

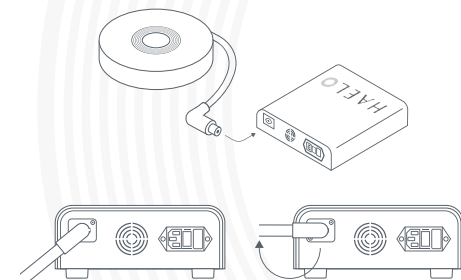

- 3 Carefully insert coil into the device, turn clockwise until the connector clicks into place. To remove, slide the silver button toward you and turn the connector counterclockwise.

- 4 Turn power switch ON.

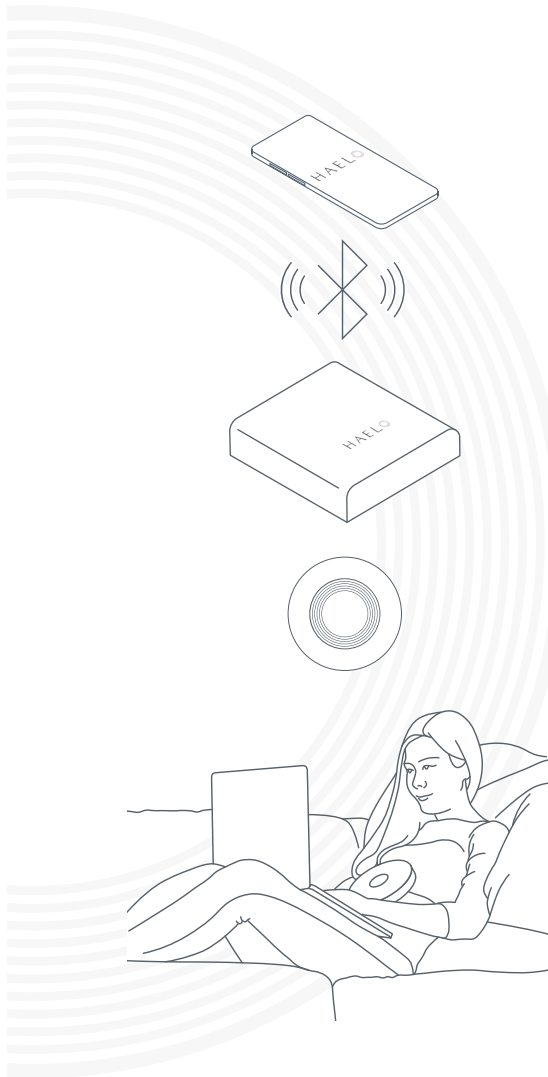

5 Download HAELO App available for iOS and Android. Your HAELO device will not work until the HAELO App is installed on your personal smart device.

6 HAELO must be powered ON to be discoverable. LEDs on the device will be illuminated. Select HAELO from the bluetooth screen.

7 Register a new account with your preferred email and password. From the App, select your desired frequency set and press PLAY!

HINT: New to PEMF? We suggest starting with the frequency set called "Harmonize." Place the coil over your stomach area, sit back and smile! Welcome to the first of many more HAELO experiences!

HAELO Frequency Sets

Ten carefully arranged frequency sets designed with you in mind!



Level Up

Powerful initiation of endorphins, serotonin, oxytocin and dopamine triggering a wave of strength, endurance and motivation to support a healthy mindset and stamina to go the distance and level up in all situations.



Prepare

Prepares the whole body and mind for exertion, activates readiness and optimized reaction time for peak performance in any challenging situation, competition or workout.



Focus

Sharpens and enhances the mind, triggers alertness, activates intuition and releases tension for restoring balance in the body and mind to outperform competition or produce an exceptional outcome.



Power

Turbo charged power not for the faint of heart! Feel like a winner and boost energy with up-leveled blood circulation, cell hydration and muscle endurance to excel physically, emotionally and mentally.



Replenish

Nourish cells with optimized hydration, circulation and oxygenation following a strenuous event, travel/jet lag or just to replenish the body. Supports the nervous system reset and helps to balance hormones.

Deep recovery of muscles, tendons, ligaments, fascia, tissue and bones after strenuous workout or competition that typically produces soreness, stiffness and pain. Supports a quick recovery for regular and consistent routines.

Recover



Spark the breath of life and align cells creating balance in the mind, organs and chakra system. Sends harmonious frequencies to enhance well being and to better support challenging and full lifestyles.

Harmonize



Calms the mind and provides stability and grounding to center the mind, body and energy field. Supports hormones and the nervous system. Targets stress reduction and reinforces peace and harmony.

Calm



Gentle soothing relief to areas that are tight or sore. Relieves physical body tension accelerating recovery. Supports optimal gut health, a sound mind and overall positive attitude to keep you on track.

Relief



Comprehensive immune defense that supports white blood cells and provides optimal protection against a wide array of nasties. Supports colon and gut health for optimum balance.

Defense



Coil Placement


The most common area for optimal coil placement for all HAELO frequency sets is the stomach or abdomen area.

Below is a guide to other suggested areas of the body for best results.

Stomach / Abdomen

Knee

Ankle

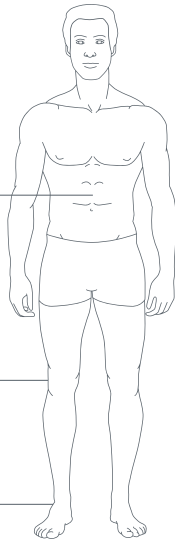

Upper Back / Shoulder

Middle Back

Lower Back



For targeted results, simply place the coil on the area of concern.



DO NOT USE THE HAELO DEVICE UNDER THESE CONDITIONS

- If you have a cardiac pacemaker, implanted defibrillator, or other implanted ferrous metallic or electronic device. Such use could cause electric shock, electrical interference or bodily harm.
- Together with a life-supporting medical electronic device such as an artificial heart, lung or respirator with ferrous material.
- Pregnant or lactating women, because the safety of PEMF devices has not been established for these conditions.
- Please check with your physician for using with children under the age of 14.
- Persons incapable of expressing their intentions or thoughts.

OUR PRODUCTS ARE NOT INTENDED FOR USE IN THE TREATMENT OR MANAGEMENT OF ANY DISEASE OR CONDITION, AND SHOULD IN NO WAY BE CONSIDERED A REPLACEMENT FOR MEDICAL ADVICE AND/OR ACTIONS TAKEN IN ORDER TO DIAGNOSE, CURE, TREAT, OR PREVENT DISEASES OF ANY NATURE AND HAVE NOT BEEN CLEARED BY THE FDA FOR ANY SUCH USE.

For complete details see Section 21

Medical Advice Disclaimer: www.haelo.com/pages/terms-of-service



Safety, Precaution & Warnings


- Avoid dropping HAELO or placing on an uneven surface or unstable location.
- Do not alter or change HAELO, as it may interfere with proper operation.
- Do not route the power cord over furniture or appliances.
- Keep power cord away from high traffic areas to avoid tripping hazards.
- Connect equipment to a grounded socket outlet and ensure the ground pin is intact. Do not defeat the ground pin.
- Do not operate HAELO with damaged or frayed power or coil cord.
- Contact HAELO Support to inquire about any product repair and replacement.
www.haelo.com/support

NOTICE Care Instructions


- Never place objects on top of the device.
- HAELO should only be cleaned with a dry cloth, avoid using excessive water or cleaning agents on the device.

HAELO Technical Specs

Front

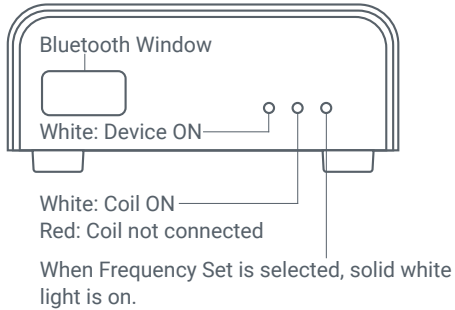

Back

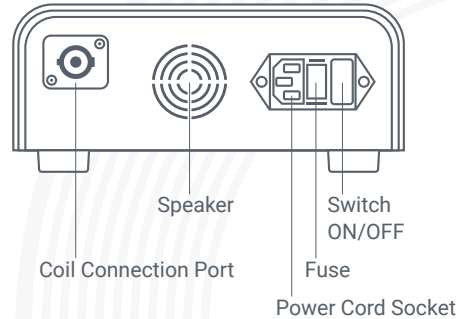

Universal AC Input	100V - 240V AC Line, 50-60Hz, 3A Fuse
Frequency Range	3Hz - 11,875Hz
Coil Rating	1.01 OHM
Minimum and Maximum Ambient Temperature Ranges	32°F to 122°F / 0°C to 50°C
<i>symphony one</i> Weight	5.5 lbs / 2.5 kgs
<i>symphony one</i> Dimensions	10.3" x 7.9" x 3.3" / 262mm x 201mm x 75mm
Coil Weight	5.5 lbs / 2.5 kgs
Coil Dimensions	8.7" Diameter x 2.6" / 221mm Diameter x 67mm
Bluetooth Low Energy (BLE)	2.4GHz Baseband

www.haelo.com/tech-specs





Limited 2-Year Warranty

HAELO warrants that for a period of two years from shipment, the device will be free from defects in material and workmanship under normal use and in accordance with appropriate handling. In the event of a defect, contact HAELO Support for instructions. www.haelo.com/support.



60-Day Money Back Guarantee

Your HAELO *symphony one* comes with a 60-day money back guarantee from date of shipment. If you're not satisfied with your purchase for any reason simply contact HAELO support for return instructions. www.haelo.com/support

www.haleo.com/legal

HAELO™ Model: *symphony one*
Input: 100-240VAC 3A 50/60Hz
S/N:



Contains: FCC ID:A8TBM64S1
IC:12246A-BM64S1

Designed and manufactured with domestic and foreign components in Silicon Valley, CA USA for haelox™ LLC. Patent Information: www.haelo.com/legal
This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Vision
Activating and optimizing every *body*.

Mission
Our mission is to improve every *body* for up-leveled performance
and optimization through the power of PEMF and
symphonic frequencies.

HAELO 

www.haelo.com/support



#howdoyouHAELO

Printed on recycled paper



20-11/061Kv1.1

